

# Brighton™

## Bipolar Coagulation Probe



- Hydrophilic Tip
- Silver spiral electrode pair design provides constant power and thermal conductivity

EXPANDING THE POSSIBILITIES™




**MERTMEDICA**  
**ENDOTEK™**

# Brighton™

## BIPOLAR COAGULATION PROBE

Merit Endotek's new Bipolar Coagulation Probe is designed with precisely the right catheter stiffness allowing for easy delivery through an endoscope's working channel while providing accurate, controlled hemostasis throughout the gastrointestinal tract. Merit Endotek offers unique probes for both single and dual plug surgical generators. Brighton Bipolar Coagulation Probes are another example of Merit Endotek's commitment to meeting the changing needs of our customers.

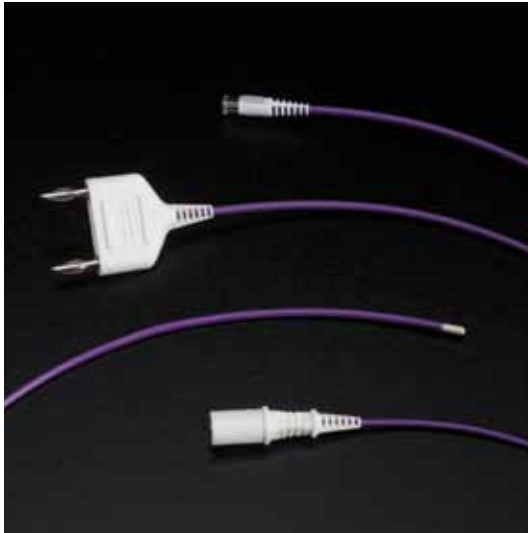


### OUR PRODUCT LINE



Fully Covered  
ALIMAXX-ES™  
Esophageal Stent



MAXXWIRE®  
Guide Wire

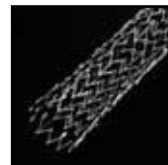

ALIMAXX-B®  
Biliary Stent



ENDOTEK BIG60™  
Inflation Device

### Brighton Bipolar Probe

Catalog Number	Guide Wire Specifications			Packaged
	DESCRIPTION	LENGTH	GENERATOR PLUG TYPE	
<b>BP-10F-DP-350</b>	10F (3.3MM)	350CM	DUAL	1 EACH
<b>BP-10F-SP-350</b>	10F (3.3MM)	350CM	SINGLE	1 EACH
<b>BP-7F-DP-350</b>	7F (2.3MM)	350CM	DUAL	1 EACH
<b>BP-7F-SP-350</b>	7F (2.3MM)	350CM	SINGLE	1 EACH

STERILE, DISPOSABLE

### BRIGHTON BIPOLAR PROBE FEATURES:

- **SILVER SPIRAL ELECTRODE PAIR DESIGN**  
provides constant power and thermal conductivity
- **HYDROPHILIC COATED TIP**  
can reduce probe adherence for more effective coagulation
- **ROUNDED DISTAL TIP**  
helps promote precise burn from various angles
- **UNIQUE CATHETER MATERIAL**  
facilitates passage through the endoscope and virtually eliminates kinking while providing optimum tamponade
- **350CM LENGTH**  
offers one size to fit most endoscopes
- **DUAL AND SINGLE PLUGS**  
to fit most bipolar generators
- **IRRIGATION THROUGH CENTRAL LUMEN**  
to effectively remove debris
- **PURPLE COLOR**  
enhances positive endoscopic identification of the tip and catheter transition



Merit Medical Systems, Inc. • 1600 West Merit Parkway • South Jordan, Utah 84095 • 1-801-253-1600 • 1-800-35-MERIT

www.merit.com

Merit Medical EUROPE, MIDDLE EAST AND AFRICA (EMEA) • Amerikalaan 42, 6199 AE Maastricht-Airport • The Netherlands • Tel: +31 43 358 82 40

Free phone for specific country:

Austria 0800 295374 • Belgium 0800 729 06 (Dutch) 0800 731 72 (French) • Denmark 80 88 00 24 • France 0800 91 60 30 • Germany 0800 182 08 71

Ireland (republic) 1800 553 163 • Luxembourg 8002 25 22 • Netherlands 0800 022 8184 • Sweden 020 79 24 45 • UK 0800 973 115 • Finland 00800 6374 8633

401985001/A ID 110410